

New Product Feature Guide



M91 TRx SL \$4995 (base model)

Super Low Single Pivot Tilt (45°)
16.75" seat to floor height
Maintains group 22 batteries
Foot platform interfaces tight to create short overall length
Standard Super Low options as per order form
Newly available options: side storage compartment; Onyx leatherette seating; SLS lighting and/or ground effects; single post Concepts back



M91 TRx Power Positioning System CG Tilt \$5395 (base model)

19.5" seat to floor height
E.A.S.E. (Easy Access Service Equipment) feature
Available with most Concepts options, including tilt, recline, extended shear reduction, Pivot Plus legs, center mount power elevating and articulating platform, Onyx seating and 17 amp/hr. battery storage for vent applications

M91 TRx Power Elevating Module \$4495

Add to your TRx CG Tilt or Recline System
20" seat to floor height,
16" x 16" – 20" x 20" seat available
7" lift = 225 lbs
6" lift = 280 lbs



SLS Single Post TRx Contour Back \$445

Sleek, new design featuring attractive, non-rehab look
Eliminates dual backcanes
Available on Super Low Tilt Systems with TRx Contour Back
Available in sizes 16" – 20" wide
5°, 10° and 15° angle adjustable back
For Onyx leatherette upgrade add \$200

Lazarus Jump Starter \$150

Helps engage the charging process of your system whenever batteries go below starting threshold of the automatic charger
Available for all systems



For more information & pricing call toll free 1-888-433-6818

Onyx Seating Package \$1140 (on-chair)

Elegant, contoured leatherette seat cover and ergonomically molded foam create an eye-catching seat and back that blurs the line between style and function
Sleek, non-rehab look
Includes automotive style headrest
Sizes from 16" x 16" – 22" x 20"




Arm Pad Storage Compartment \$175 each

Ingenious armrest with handy flip up storage box
Complete with slide out cup holder
Available only on TRx Power Positioning Systems

Armrest Storage Compartment \$125 each

First ever, easy to operate storage compartment
Sleek, modern, automotive look
Special release mechanism for clients with limited dexterity
Ideal for storing personal items
Available only on TRx Tilt and Tilt/Elevating Systems

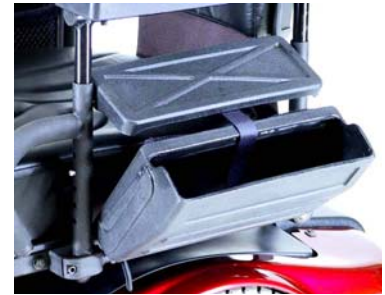

Latitude \$7995 (base model)

An amazing, versatile system designed so clients can easily maneuver themselves from floor, to seating level and beyond
12" x 12" – 15" x 15" seat
175 lb. weight capacity
Height, width and depth adjustable with 5" of seat-to-floor height adjustability
Available only on Storm short power bases

Anterior Tilt Module \$3995

Up to 30° anterior tilt
Maintains 55° of posterior tilt
Available on power elevating CG tilt systems
Available on 626, 646, Storm and select Quantum power bases
Aids in forward standing transfers
Power recline module can be added



For more information & pricing call toll free 1-888-433-6818

Jazzy SLS Ground Effects \$295

Stylish, rear spoiler upgrade package
Visually appealing
Front ATS wheel covers



SLS Lighting Effects \$350

Exclusive package featuring amazing, eye-catching neon lights
Available on any tilt system 16" x 16" and greater
(shown on Jazzy SLS System)



Pivot Plus Power Elevating Foot Platform \$2695

No need for articulation!
Unique, durable design allows entire platform to be adjusted in/out and up/down to align with client's knee pivot (co-axial symmetry)
Removable, weight bearing construction
Convenient flip up footplate
Adjustable for most leg lengths
Available on 16" x 16" or greater

Center Mount Power Articulating & Elevating Platform \$2995

Extends as it elevates and articulates
Center mount has a shorter swing than standard front rigging
Flip up footplate
Excellent for side transfers
Articulates 5.5"
Available on 18" standard or greater



Manual Pivot Plus Legs \$795 (on-chair)

No need for articulation!
Unique design allows legrest to be adjusted in/out and up/down to align with client's knee pivot (co-axial symmetry)
Heavy duty rubber dampening system and footplate rollers assist in opening doors
Footplates designed for minimal system length while maintaining caster clearance
Adjustable for most leg lengths

1121 TRx SL \$4995 (base model)

Featuring a 16.75" seat to floor height
Uses the same great system options available
on the 1122 Super Low
16" x 16" – 20" x 20" seat size; 250 lb. maximum



TRx Contour Arm Troughs \$85 each Hand pad available \$55 each

Highly attractive, ergonomic arm trough pads
Keeps clients arm steady in any position while tilting
Avoids digging into elbows while reclining
Multi position cutouts for versatile arm
position flexibility

New Flip Back Tilt & Reclining Arms

Virtually eliminates any play
New flip back auto stop feature
Eliminates the risk of injuring the joystick as
it rests backwards
Standard on all systems effective January 1, 2003



M44 'R' Scanner Electronics Box \$495

Amazing scanner electronics shows clients
what mode their chair is currently in
Designed to operate a multi function system with a
single momentary control switch
Ideal for quadriplegics or those with limited
hand dexterity

'Precline' \$1500

Snugly, securely and comfortably
'sandwiches' clients in a lean, forward position
Up to 30° forward recline available
Range of recline is decreased by degrees of
precline, limited to 173° total
0° – 8° available at no charge



U.S.
700 Ensminger Rd., Suite 112
Tonawanda, NY 14150
Toll Free Tel: 1.888.433.6818
Toll Free Fax: 1.888.433.6834

CANADA
84 Citation Dr.
Concord, ON L4K 3C1
Tel: 905.695.0134
Fax: 905.695.0138