

ADDITIONAL MEASUREMENTS FOR CLIENTS WITH LYMPHEDEMA


Dealer Name: _____ Dealer Account #: _____

PO #: _____ Tag: _____


Client's Name: _____ Client's current weight: _____

Client's Weight History: _____

Left Lateral View:

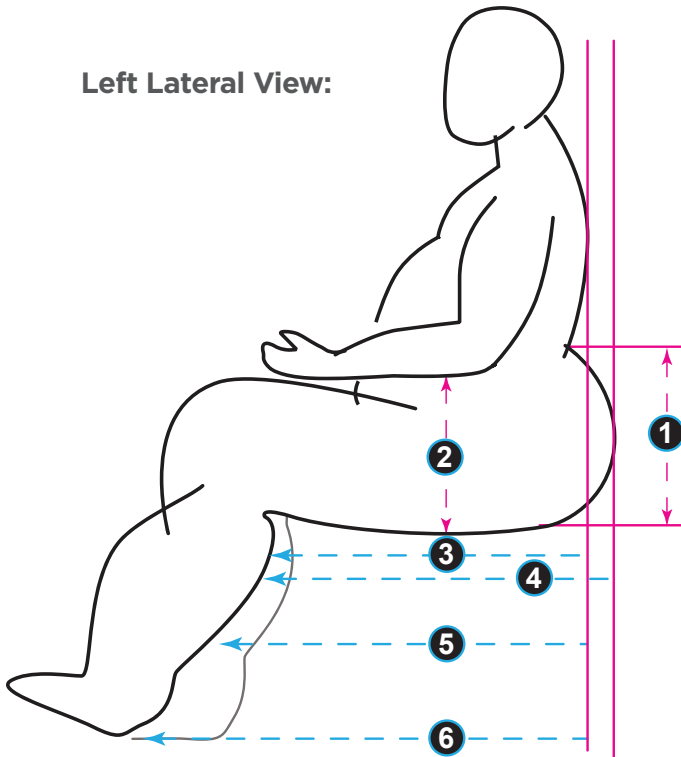

Right Lateral View:

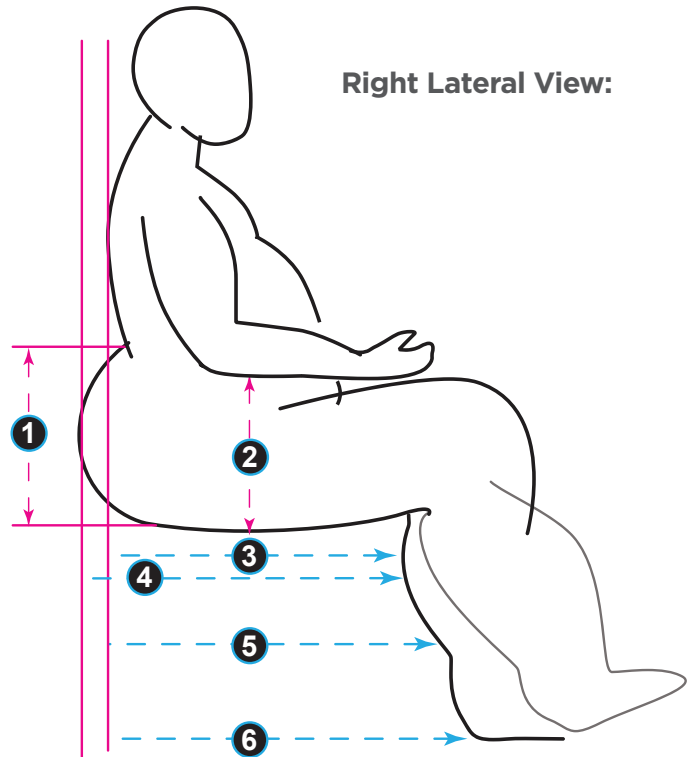


1. **Buttocks height:** vertical distance from the top of the sitting surface to the top of gluteal tissue
Left side: _____ Right side: _____
2. **Elbow height:** vertical distance from the top of the sitting surface to the inferior aspect of the elbow in the sitting position
Left side: _____ Right side: _____
3. **Sitting depth:** distance from the most posterior aspect of the trunk to tissue behind the knee (for functional seat depth)
Left side: _____ Right side: _____

4. **Buttocks/thigh depth:** distance from most posterior aspect of the buttocks to the most posterior aspect of tissue behind the knee
Left side: _____ Right side: _____
5. **Calf depth:** distance from the most posterior aspect of the posterior aspect of tissue behind the knee to the most posterior aspect of the calf tissue
Left side: _____ Right side: _____
6. **Lower extremity depth:** distance from the most posterior aspect of tissue behind the knee to the posterior heel
Left side: _____ Right side: _____

Right Side:

Left Side:

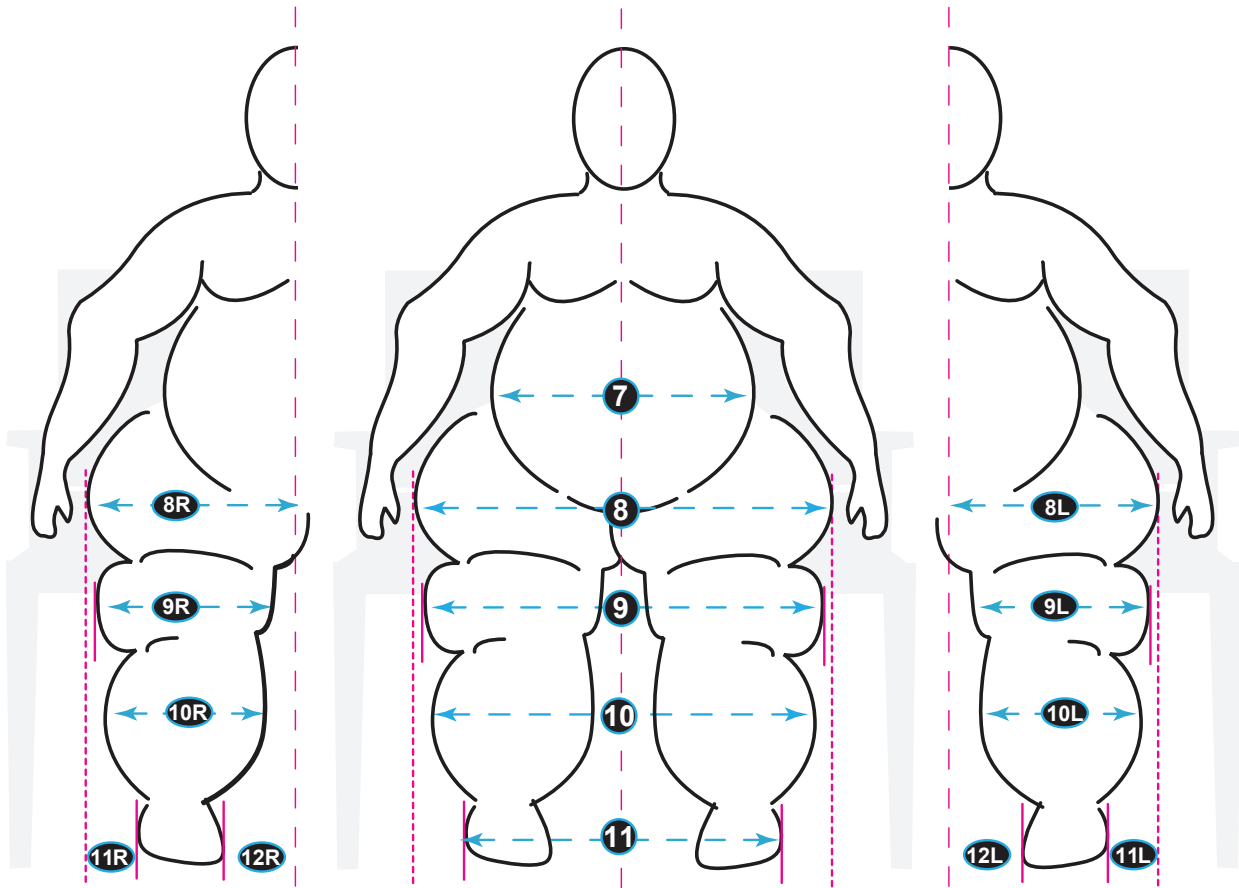

7. **Trunk width:** overall width of the trunk, at lower ribs or widest point _____

8. **External thigh width:** width at widest point (used to determine seat width) _____

Left side: _____ Right side: _____

9. **External knee width:** width at widest point lateral left knee tissue to lateral right knee tissue _____

Left side: _____ Right side: _____

10. **External lower leg width:** distance between the outside borders of the lower legs at their widest point _____

Left side: _____ Right side: _____

11. **External foot width:** distance between the outside borders of the feet _____

12. **Lateral thigh to lateral foot:** Horizontal distance between most lateral aspect of the thigh to the lateral aspect of the foot

Left side: _____ Right side: _____

13. **Medial foot to centerline:** Distance from the medial aspect of the foot to centerline

Left side: _____ Right side: _____